

# Safety blood collection sets

## NEW BD Vacutainer® UltraTouch™ Push Button blood collection set

Further to the features of the BD Vacutainer® Push Button Blood Collection, the BD Vacutainer® UltraTouch collection set employs patented PentaPoint™ Comfort 5-bevel needle technology. Studies have shown that this design helps reduce the chance of a painful injection by creating a flatter, thinner surface to help penetrate the skin with significantly greater ease.<sup>1</sup> This new safety device has been shown to reduce penetration forces by up to 32 percent when compared to another leading blood collection set.<sup>2</sup>

In addition, its exclusive BD RightGauge™, Ultra-thin wall technology allows a better blood flow due to the needle's larger inner diameter. Therefore, clinicians can select a smaller gauge needle without sacrificing sample quality. This technology can also improve tube fill time by up to 50% when using the same gauge (23G standard needle vs. 23G UltraTouch™).<sup>3</sup>

UltraTouch™ Push Button brings healthcare workers greater confidence and ability in knowing they can use the needle gauge that is most appropriate for their patients with improved efficiency.


### Further information

BD recommends for all needles without an integrated holder that a BD Vacutainer® One Use holder is used.

Clinical and technical information is available on request.



**NEW**

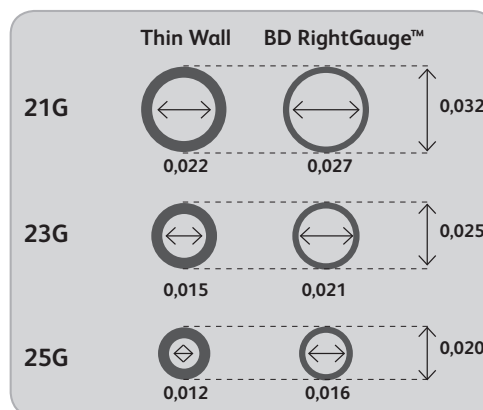

All dimensions shown are nominal values and in inches.

## BD Vacutainer® UltraTouch™ Push Button blood collection sets

Cat. no.	Size	Needle length	Length of tubing	With Luer adapter	Colour code	Box/case
367393	21G (0.8mm)	19mm	178mm	Yes		50/200
367365	21G (0.8mm)	19mm	305mm	Yes		50/200
367392	23G (0.6mm)	19mm	178mm	Yes		50/200
367364	23G (0.6mm)	19mm	305mm	Yes		50/200
367391	25G (0.5mm)	19mm	178mm	Yes		50/200
367363	25G (0.5mm)	19mm	305mm	Yes		50/200

1. Hirsch L.J., et al. Journal of Diabetes Science and Technology. 2012, 6(2):328-35.  
 2. BD White Paper VS9248: BD Vacutainer® Push Button Blood Collection Set Penetration Force, 2016.  
 3. BD White Paper VS9249: BD Vacutainer® Push Button Blood Collection Set Tube Fill Time, 2016.